

WHERE ESTHETICS MEETS FUNCTION

ARTICULATOR SYSTEMS

2026 PRODUCT COLLECTION



PANADENT CORPORATION

TO ORDER

EMAIL: INFO@PANADENT.COM

INSTAGRAM: [@PANADENTCORP](https://www.instagram.com/PANADENTCORP)


PHONE: (800) 368-9777

BEST SELLING ARTICULATORS & PACKAGES



1610 Magnetic Model PCH Articulator

The model PCH articulator is a precision articulator which allows the transfer of dental casts from one PCH or PSH articulator to another. The model PCH has a curved pin calibrated in degrees. The 1610 AR Comes with 20 Magnetic Mounting Plates & Hex Wrench.

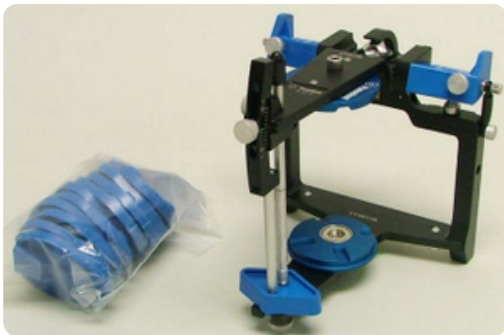

1650 Magnetic Model PCF Articulator

The model PCF articulator is a precision articulator featuring fixed 1.5mm analogs set at 20°, an installed magnetic system and a curved incisal pin calibrated in degrees. The 1650 AR Comes with 10 Magnetic Mounting Plates.


1210 Basic Module with PCH Articulator



The Basic Module includes all of the essential items listed below needed to mount casts when using the Pana-Mount™ Face-bow. The molded plastic case holds the articulator, face-bow, and accessories plus room for mounted casts.

Includes: 1610 Magnetic Model PCH Articulator, 2094 Articulator Support Legs (set), 4000 Pana-Mount™ Face-bow, 4070 Bio-Esthetic™ Level Gauge, 8210 Molded Plastic Case, 2280 Mandibular Mounting Stand, 7172 Bite-Tray™ Registration Plates (20), and 4410 Bite-Fork Stabilizers (60).


DIGITAL & ANALOG ACCESSORIES

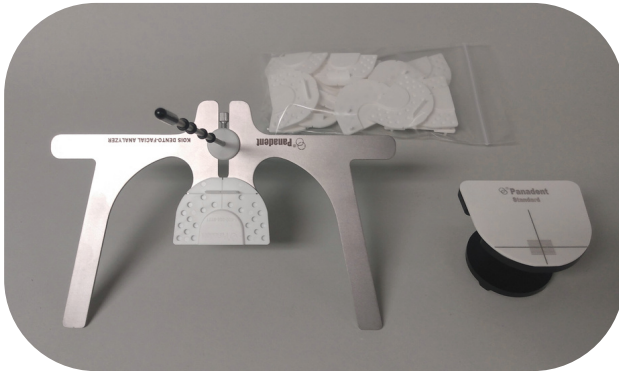

4300 KOIS ANALYZER SYSTEM

The Kois Dento-Facial Analyzer registers and transfers the patient's occlusal plane in three planes of space to the articulator. Based on Dr. John Kois 100mm axis-incisal distance research. Includes: Analyzer, Platform, and 20 Index Trays.

4305 Kois Dento-Facial Analyzer and 4360 Kois Platform also available separately.

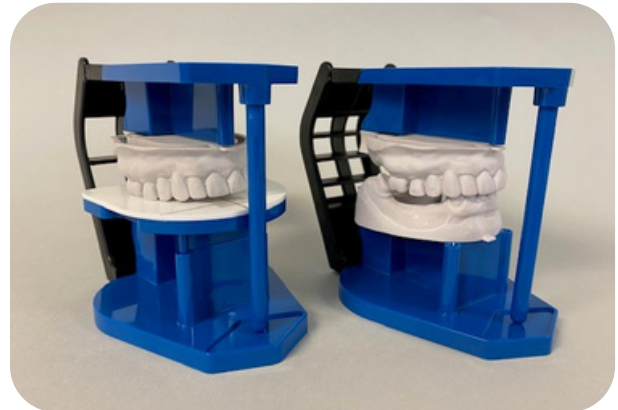

1500 PAL 2.0 MINI 2-IN-1 ARTICULATOR

Panadent PAL 2.0 Articulator System is a reusable versatile 2in1 plastic articulator that has the same axis as the full size Panadent articulator designed as an independent system. It can mount traditional models with plaster and 3D printed models without plaster. Kit includes blue plastic upper/lower connectors, plastic platform, and PAL.



4000 PANA-MOUNT FACE-BOW


The Pana-mount Facebow System orients the upper cast on the articulator in the same spatial relationship as the patient's maxilla to their TMJ. Includes: Face-Bow, two complete Bite-Fork Assemblies, Mounting Fixture, Standard Nasion Relator, Screws, and Wrench



4380 KOIS DIGITAL TRANSFER ADAPTER

The Kois Digital Transfer Adapter is used to mount 3D printed models without plaster on all Magnetic Panadent PCH, PCF, and PSH articulators using Exocad, Blender, and Meshmixer software. A great way to simplify your workflow!

ESTHETIC & OCCLUSION TOOLS

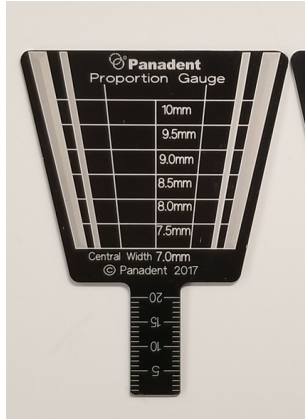
9375


ESTHETIC TOOTH PROPORTION GAUGE

Designed to provide "Esthetic Zones" for tooth widths and lengths. The tooth width proportion "Esthetic Zones" range from .618 to .75 which includes Golden Proportion as well as RED (Reoccurring Esthetic Dental) Proportion. The tooth length proportion "Esthetic Zone" chart ranges from 75% to 85% width to height ratios. Made of aluminum & is autoclavable.


9391

MULTI-MEASURING RULER

The MMR is used for measuring tooth lengths, tooth widths, vertical overbite, horizontal overjet, as well as having opening and lateral range of motion scales. Made of durable aluminum and is autoclavable. The perfect mm ruler for your practice.

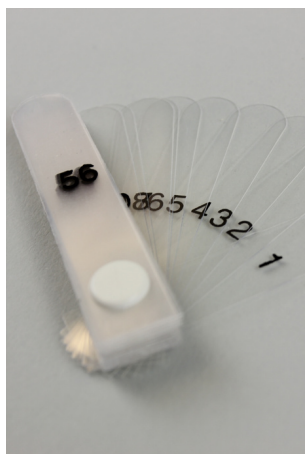
9390

LEAF GAUGE (NUMBERED)

This leaf gauge is made of polyester strips and are numbered (1-56) in 1/10th mm increments for convenience. They are used as an anterior deprogrammer for taking centric records and can also be used in equilibration or coronaplasty procedures.